

# HI-EXPECTATIONS MULTIFOCAL

Service to the Eyecare Professions since 1902

## HI-EXPECTATIONS MULTIFOCAL

**Hi-Expectations** is Winchester Optical's unique tri-aspherical multifocal gas permeable lens design. **Hi-Expectations** incorporates a back bi-aspherical design, along with a front spherical/aspherical geometry for high add requirements.




**Hi-Expectations** uses our original window design allowing complete control over both the distance and near zones. The larger the distance window the stronger the concentration of distance power, conversely the smaller the distance window the stronger the concentration of near power.



All **Hi-Expectations** multifocals are made in Boston ES®, advanced Aercor® chemistry which transfers excellent oxygen permeability without compromising surface or structural character. The high oxygen permeability of Boston ES assures your patients' corneal health and unsurpassed comfort.

HI-EXPECTATIONS PARAMETERS	
<b>POWER</b>	+20.00 to -20.00
<b>BASE CURVES</b>	7.00 mm to 8.40 mm in .10mm steps
<b>DIAMETER</b>	9.6
<b>ADD POWERS</b>	up to +3.75

Hi EXPECTATIONS TRIAL SET SPECIAL OFFER	
11 Lens Diagnostic Trial Set	

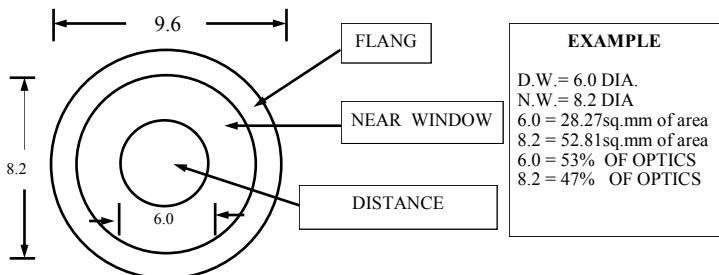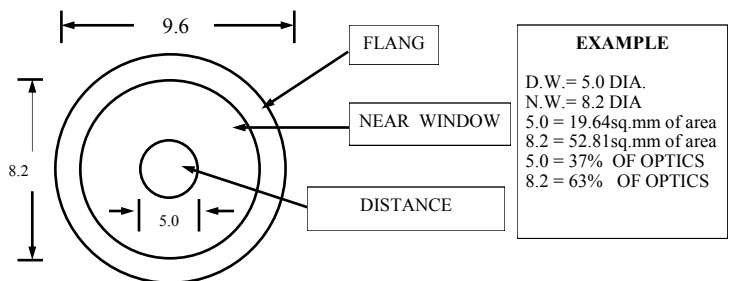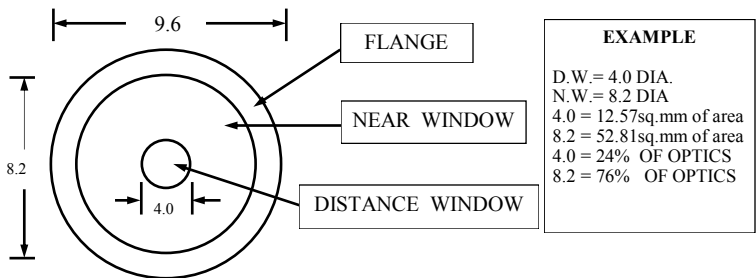

**OUR SERVICE GUARANTEE ON ALL HI-EXPECTATIONS  
 --LENSES ORDERED BY 11:00 AM WILL BE SHIPPED THE SAME DAY--**